


Bode RapidSwab® Instructions for Use-Sample Collection


FOR HUMAN/FORENSIC IDENTIFICATION ONLY—NOT FOR CLINICAL OR DIAGNOSTIC USE

Note: Collection personnel should wear proper personal protective equipment at all times during collection of DNA evidence.

COLLECTOR COMPONENTS

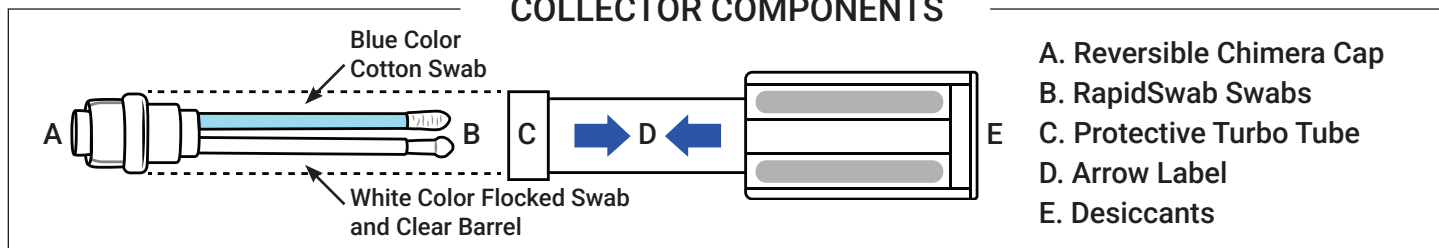


Step 1

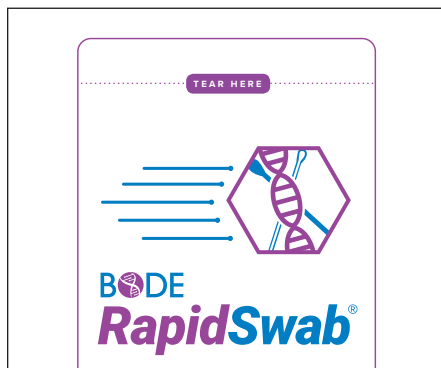

Tear open the foil wrapper at the notches and remove the collector holding the protective Turbo Tube. If needed, remove a barcode from the back of the pouch and apply it to the Turbo Tube.


The remaining barcodes can be used as needed for labeling Transport Pouches, envelopes, or other documentation.

Step 2

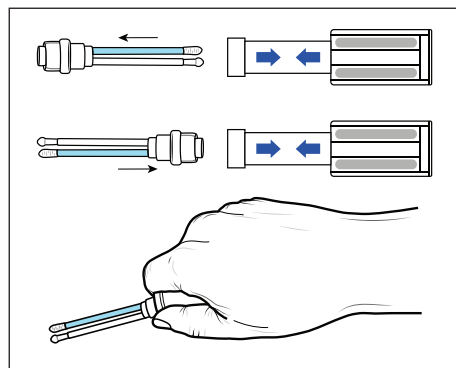


Carefully remove the reversible Chimera Cap with attached swabs from the Turbo Tube and insert into the Turbo Tube as shown above. Do not touch the swab heads when handling the device.

Step 3

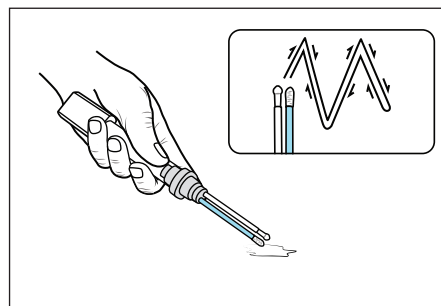


- For Dry Stain Collection:** Add a few drops of water (2-3) to one side of the cotton swab (blue swab stick) and a few drops (1-2) to one side of the flocked swab. Swab over the biological stain/evidence item in an "M" motion (representation above) to simultaneously collect with both swabs (wet side) before rotating and swabbing again with the dry side.
- For a Wet Stain Collection:** Swab the biological stain/evidence item in an "M" motion (representation above) to simultaneously collect with both swabs. Rotate the collector and swabs as needed to collect as much of the stain as possible.

Step 4

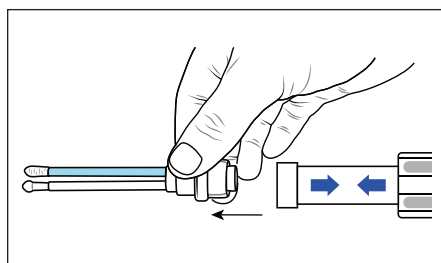


Remove the reversible Chimera Cap from the Turbo Tube by holding the cap firmly and twisting the Turbo Tube/Cap to separate. Be careful not to touch the swabs during this step.

Step 5

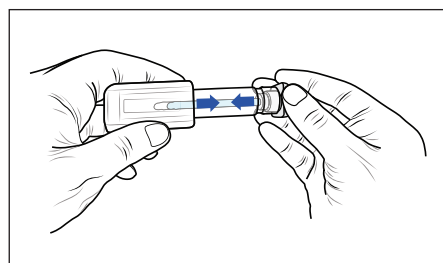

Insert the swabs back into the Turbo Tube aligning the blue swab stick in between the arrow label. When inserted correctly, you can visually see the blue swab stick as pictured above.

Step 6

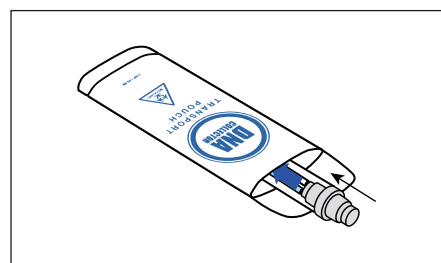

Insert the RapidSwab Collector into evidence packaging such as a Bode Transport Pouch. Seal and label per your agency's requirements. The foil wrapper and remaining unused barcodes can be discarded.




Bode RapidSwab® Instructions for Use-Sample Analysis


FOR HUMAN/FORENSIC IDENTIFICATION ONLY—NOT FOR CLINICAL OR DIAGNOSTIC USE

Note: Collection personnel should wear proper personal protective equipment at all times during analysis of DNA evidence.

COLLECTOR COMPONENTS

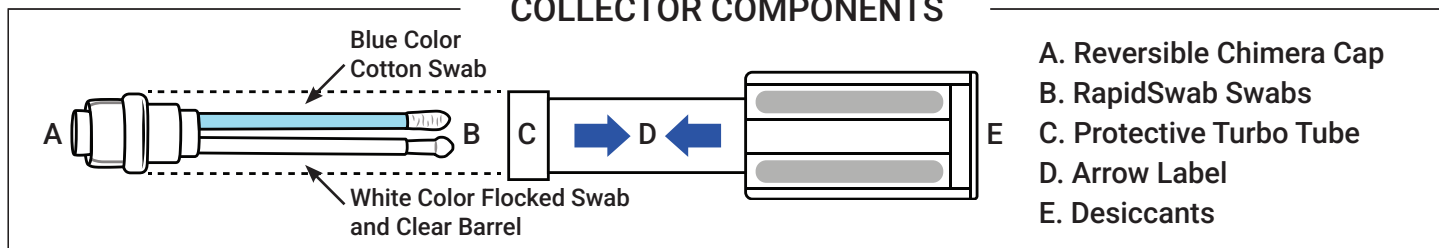


Step 1

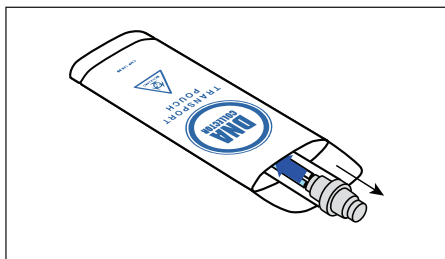


Open the Transport Pouch or packaging containing the collected RapidSwab Sample.

Step 2

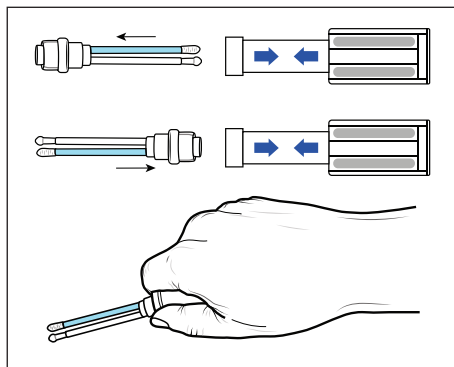


Carefully remove the reversible Chimera Cap with attached swabs from the Turbo Tube and insert into the Turbo Tube as shown above. Do not touch the swab heads when handling the device.

Step 3

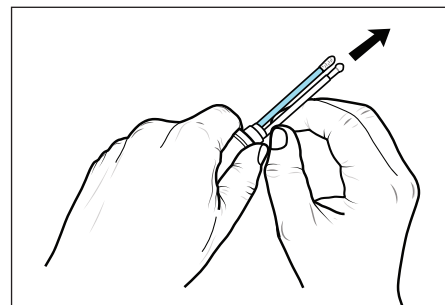


Holding the Chimera Cap in one hand, firmly grasp the clear barrel and flocked swab. Gently remove the Flocked Swab and barrel from the Chimera Cap by pulling up and away from the Chimera Cap.

Step 4

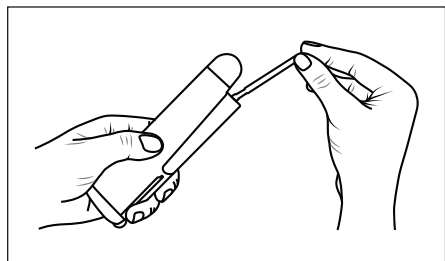


Insert the flocked swab into the rapid cartridge, the barrel can be discarded.

Step 5

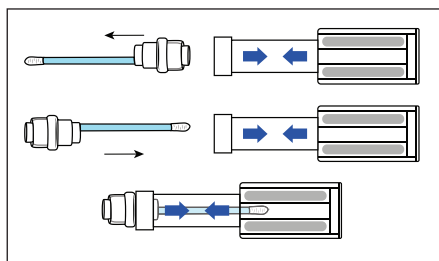

Remove the reversible Chimera Cap from the Turbo Tube by holding the cap firmly and twisting the Turbo Tube/Cap to separate. Insert the Chimera Cap back into the Turbo Tube aligning the blue swab stick in between the arrow label. When inserted correctly, you can visually see the blue swab stick as pictured above.

Step 6

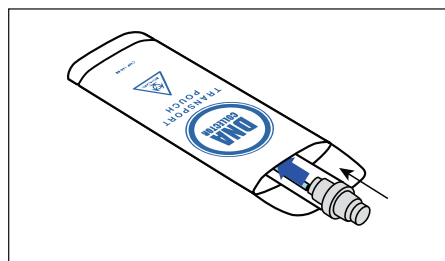

Repackage the RapidSwab containing only the cotton swab into the original Transport Pouch/Evidence Packaging.

